

CHAPTERS

Tables/
Breadboards

Mechanics

Optomechanic
Devices

Kits

Lab Supplies

SECTIONS

Mini Series

Ø1/2" Post
AssembliesØ1" Post
AssembliesØ1.5" Post
AssembliesMounting/Angle
Brackets

Lab Platforms

Instrument Shelves

Lens Tubes

Cage Systems

Optical Rails

Mirror Mounts

Kinematic Mounts

Fixed Mounts

Translation Mounts

Rotation Mounts

Filter Mounts

V-Mounts


Iris/Apertures

Adapters

SM2-Compatible Quick-Release Lens Tube Couplers

Also Available in
Ø1/2", Ø1.00",
and Ø3.00"

SM2Q

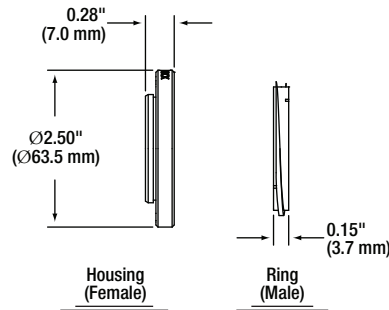

Please refer to our website for complete models and drawings.

Easily Swap
Components

- One Full Turn Connect and Disconnect Adapters
- Lockable Design
- Easily Swap SM2-Threaded Components

The SM2Q Adapter Set provides a means to quickly connect and disconnect two similarly SM2-threaded components such as lens tubes, mounted optics, and detectors. Thread the male adapter onto an externally SM2-threaded component, and thread the female adapter onto an internally SM2-threaded component. A witness line on each adapter allows them to be aligned easily.


ITEM #	\$	£	€	RMB	DESCRIPTION
SM2Q	\$ 48.00	£ 34.56	€ 41,76	¥ 382.56	SM2 Quick-Release Adapter Set

SM2-Compatible, Lever-Actuated Iris Diaphragm

Also Available in
Ø1/2",
Ø1.00", and Ø3.00"

SM2D25

Our standard lever-actuated iris diaphragm is continuously variable and compatible with our Ø2.00" SM2-threaded lens tubes. It has been designed for reliability over thousands of cycles.

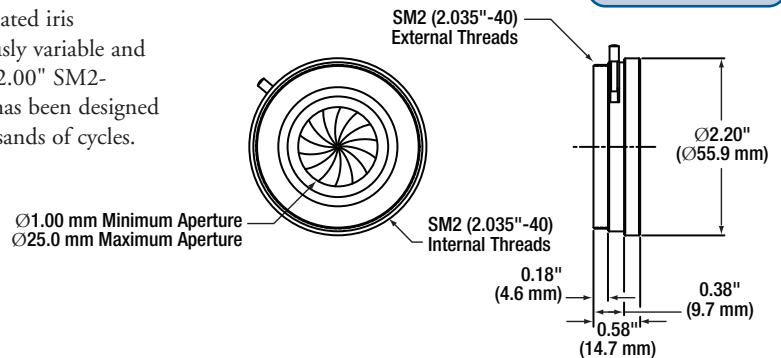See Pages 323 - 327 for Our
Extensive Line of Irides


ITEM #	\$	£	€	RMB	APERTURE RANGE	DESCRIPTION
SM2D25	\$ 84.70	£ 60.98	€ 73,69	¥ 675.06	Ø1.0 mm to Ø25.0 mm	SM2 Lever-Activated Iris Diaphragm

SM2-Compatible, Ring-Activated Iris Diaphragm

Also Available in
Ø1/2" and
Ø1.00"

SM2D25D

- Ring-Activated Iris Diaphragm Ideal for Cage System Applications
- Design Eliminates Actuating Arm, Significantly Reducing Clearance Requirements

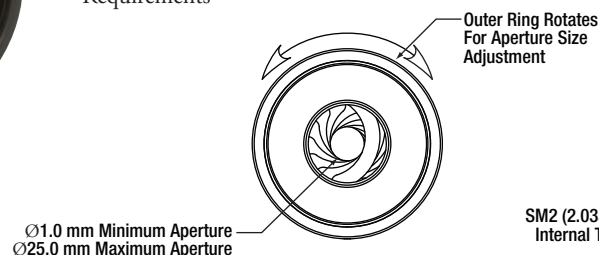

Mechanical
Drawings Available on the
WEB

ITEM #	\$	£	€	RMB	APERTURE RANGE	DESCRIPTION
SM2D25D	\$ 75.50	£ 54.36	€ 65,69	¥ 601.74	Ø1.0 mm to Ø25.0 mm	SM2 Ring-Activated Iris Diaphragm