

## Microscope Adapters with Internal SM Series Threads


### Features

- Integrate SM1, SM2, and 30 mm Cage System Products with Microscopes
- Adapters for the Following Microscopes:
  - Olympus IX and BX Microscopes
  - Nikon Eclipse and Eclipse Ti Mount Microscopes
  - Leica DMI Microscopes
  - Zeiss Axioskop Microscopes

These adapters allow Olympus BX and IX, Leica DMI, Nikon Eclipse, Nikon Eclipse Ti, or Zeiss Axioskop microscopes to be mounted with SM Series lens tubes and integrated with our cage assemblies.

The adapters mate to the reflected light lamphouse port of current generation microscope models. Any existing lamphouse needs to be removed from the microscope prior to using the adapter.

**NEW**  
versions

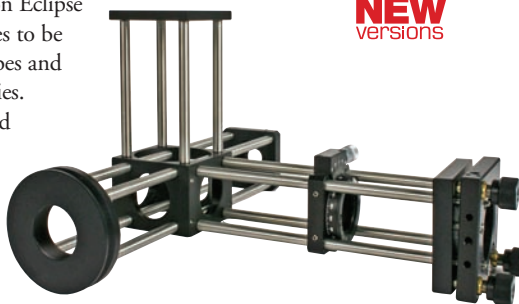

COMPATABILITY	OLYMPUS BX & IX ITEM #	LEICA DMI ITEM #	ZEISS AXIOSKOP ITEM #	NIKON ECLIPSE ITEM #	NIKON ECLIPSE TI ITEM #
Internal SM1 (1.035" -40) Thread and 30 mm Cage	SM1A14	SM1A21	SM1A23	SM1A22	SM1A26*



COMPATABILITY	OLYMPUS BX & IX ITEM #	LEICA DMI ITEM #	ZEISS AXIOSKOP ITEM #	NIKON ECLIPSE ITEM #	NIKON ECLIPSE TI ITEM #
Internal SM2 (2.035" -40) Thread	SM2A7	SM2A8	SM2A10	SM2A9	SM2A18*

### Adapters with Internal SM1 Threads, 30 mm Cage System Compatible

ITEM #	\$	£	€	RMB	DESCRIPTION
SM1A14	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Olympus BX and IX to Internal SM1 Thread, 30 mm Cage System Adapter
SM1A21	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Leica DMI to Internal SM1 Thread, 30 mm Cage System Adapter
SM1A23	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Zeiss Axioskop to Internal SM1 Thread, 30 mm Cage System Adapter
SM1A22	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Nikon Eclipse to Internal SM1 Thread, 30 mm Cage System Adapter
SM1A26*	\$ 110.11	£ 79.28	€ 95,80	¥ 877.58	Nikon Eclipse Ti to Internal SM1 Thread, 30 mm Cage System Adapter

\*These Nikon Eclipse Ti bayonet adapters are the same as the Nikon Eclipse adapters except that they incorporate an additional spring. Whether or not you need an adapter with a spring will depend on the compatibility requirements of the light port on your Nikon microscope.

### Adapters with Internal SM2 Threads

ITEM #	\$	£	€	RMB	DESCRIPTION
SM2A7	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Olympus BX and IX to Internal SM2 Thread Adapter
SM2A8	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Leica DMI to Internal SM2 Thread Adapter
SM2A10	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Zeiss Axioskop to Internal SM2 Thread Adapter
SM2A9	\$ 88.00	£ 63.36	€ 76,56	¥ 701.36	Nikon Eclipse to Internal SM2 Thread Adapter
SM2A18*	\$ 110.11	£ 79.28	€ 95,80	¥ 877.58	Nikon Eclipse Ti to Internal SM2 Thread Adapter

\*These Nikon Eclipse Ti bayonet adapters are the same as the Nikon Eclipse adapters except that they incorporate an additional spring. Whether or not you need an adapter with a spring will depend on the compatibility requirements of the light port on your Nikon microscope.