

Keyboard Holder

- Includes Wrist Pad
- Attaches to ScienceDesk Frame


The frame-mounted Keyboard Holder is suitable for most keyboards and mouse pads. It fixes to the side rails and provides 360° of horizontal positioning and vertical adjustment.

Keyboard Holder

ITEM#	\$	£	€	RMB	DESCRIPTION
PSY132	\$ 399.00	£ 276.60	€ 354.30	¥ 3,369.20	Frame-Mounted Keyboard Support



Flat Screen Support Bracket

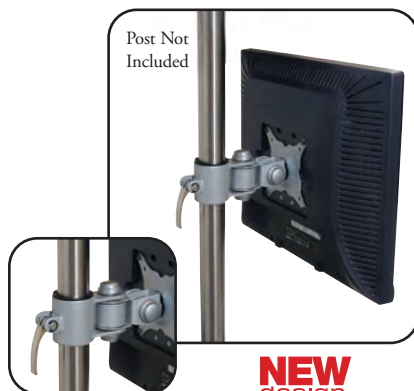

- Compatible with 75 mm x 75 mm and 100 mm x 100 mm VESA Patterns
- Mounts to a Ø1.5" Post

The flat screen support bracket is compatible with 75 mm x 75 mm and 100 mm x 100 mm VESA mounting patterns. These mounting patterns are typically used on small to medium LCD displays. A Ø1.5" post must be used with this mount and is not included. An acceptable post is our PSY160 (below); alternatively, a post from an Overhead Shelf or Instrument Shelf may also be used.



Flat Screen Support Bracket

ITEM#	\$	£	€	RMB	DESCRIPTION
PSY120	\$ 105.00	£ 72.80	€ 93.30	¥ 886.70	Flat Screen Support Bracket

Post Assembly



- Total Length: 24" (600 mm)
- Usable Length: 23.5" (588 mm)
- Ideal for Mounting Ø1.5" Post Accessories
- Mounting Hardware Included

These posts are suitable for mounting auxiliary equipment. They can be ergonomically positioned and mounted on the surrounding frame of the ScienceDesk. The post assembly is offered with a 1.5" (38 mm) outer diameter and a usable length of 23.5" (588 mm).

Post Assembly

ITEM#	\$	£	€	RMB	SIZE
PSY160	\$ 65.00	£ 45.10	€ 57.80	¥ 548.90	Ø1.5" x 24" (Ø38 mm x 600 mm)


Post Mount Adapter

- Mount a Ø1.5" ScienceDesk Post (PSY160) to a Breadboard
- Metric and Imperial Compatible
- Use ScienceDesk Accessories on an Optical Table

The TA1 post adapter allows a Ø1.5" ScienceDesk post (PSY160) to attach to a breadboard or optical table. The adapter has through holes spaced 50 mm x 50 mm (2" x 2") to attach to standard metric (M6 x 1.0) or imperial (1/4"-20) breadboards. With this adapter, the ScienceDesk Accessory System can compliment nearly any standard breadboard, providing increased usable work space. The ScienceDesk posts attach to the adapter via an M16 counterbore. The adapter includes necessary metric and imperial mounting hardware.

Please note that this adapter is not designed for use with two post shelves such as the Overhead Shelves.

**NEW
product**

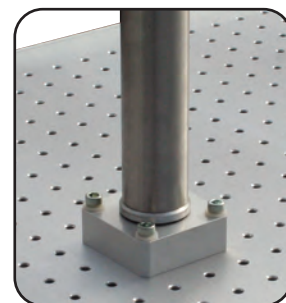

TA1 shown with PSY160
(Sold Separately)

Post Mount Adapter

ITEM#	\$	£	€	RMB	SIZE
TA1	\$ 27.50	£ 19.10	€ 24.50	¥ 232.30	Post Mount Adapter